

Miscellaneous Products

Barraquer Wire Speculum

8207

Barraquer Wire Speculum
Solid blades, reusable

Adult (matte finish)
Blade Width: 14 mm
Blade Spread: 19 mm



Phoropter Near Vision Rod

93692

Holds Phoropter clip with near vision card. Fits most standard phoropters. Length: 26.5"




ROD ONLY #93692

VU-Filterz® Occlusion Sheet

1 sheet/pkg, size: 9.5" x 4.1" (240 x 105 mm)
Up to 5 filters can be cut out of one sheet depending on the size of the spectacle lens.
(Compare to Bangerter Foils)

VU-Filterz is an effective substitute to traditional occlusion devices for both children and adults. VU-Filterz is an extremely thin polymer available in graduated densities that range from a very subtle reduction of visual acuity to a grey filter that offers total occlusion. The filters are cut to fit a spectacle lens by the eye care professional and are easily applied in minutes. For adults, VU-Filterz are a cosmetically acceptable means of eliminating disturbing diplopia at the higher densities.

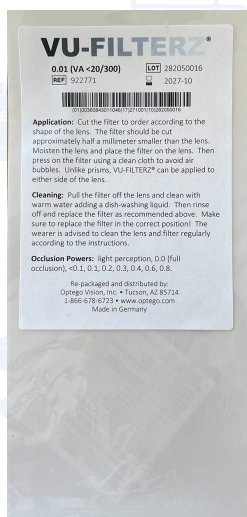


1 VU-Filterz Sheet/pkg up to 5 filters per sheet

0.8 (VA 20/25)	922111
0.6 (VA 20/30)	922221
0.4 (VA 20/50)	922331
0.3 (VA 20/70)	922441
0.2 (VA 20/100)	922551
0.1 VA 20/200)	922661
0.01 (VA <20/300)	922771
Light perception	922881
0.0 (VA none)	922991

Tonometer Tip Cleaning Station

95385

Tonometer Tip Cleaning Station for GOLDMANN tonometer prisms with screw cap and drying ports
Holds 6 tonometer tips Diameter: 3-1/8" (80 mm)

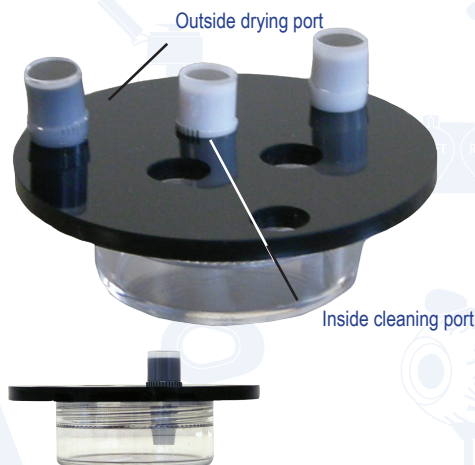

Outside drying port

Inside cleaning port

The tonometer tip cleaning station #95385 has a screw cap and six ports to hold the tonometer tips. The top unscrews so that cleaning/disinfection solution may be added to the base. The three inside cleaning ports hold the used tonometer tips in the cleaning solution. Once cleaned and disinfected, the tips are moved to one of the three outside drying ports.

Chin Rest Paper Pins

Stainless steel 2/pkg
Universal, will fit most slit lamps

98425



Plastic 2/pkg
Universal, will fit most slit lamps

98428



Alcohol Dispenser

Liquid well pump
7.6 fl.oz. (225 ml)

9949

Excellent for dispensing isopropyl alcohol. The clear bottle allows for fluid level visibility. Unique push design, preventing cross contamination.



OPTEGO VISION USA
4525 South Coach Dr
Tucson, AZ 85714
Phone 866-678-6723
Fax 520-209-1586
email@optego.com
optegousa.com

OPTEGO VISION ASIA
60 Paya Lebar Rd #11-33
Paya Lebar Square
Singapore 409051
Phone +65 6742 5940
asia@optego.com
optego.com.sg