
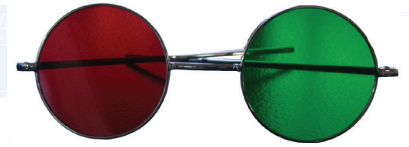


Worthmore® II Four Dot Test

94195

Red-Green Glasses - reversible

9397

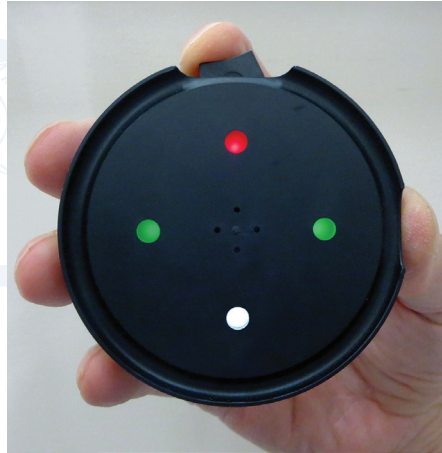

Red-Green Reversible Glasses
Chrome frame

Red-Green Glasses - adult

1 pair
2 pairs

9374
93748



Plastic, 1 pair, adult size
Red lens on right side - green lens on left



Plastic, 2 pairs, adult size
1 pair red lens on right - green lens on left
1 pair red lens on left - green lens on right

Red-Green Glasses - pediatric size

9375



Plastic, 1 pair, pediatric size
Frame color may vary

Worthmore® II Four Dot Test with a pair of red-green glasses (#9397). Lightweight, palm-sized.

This ergonomically designed handheld device features 4 peripheral porthole LEDs (Light-Emitting Diodes) and 4 central porthole LEDs with a constant brightness. The peripheral porthole lights test for peripheral suppression; the central porthole lights test for foveal fusion/suppression.

A momentary rocker switch allows for peripheral or central testing.

The Worthmore II Four Dot Test is powered by two replaceable coin batteries (included)
Size: Diameter 2.5" / Thickness: 0.5"

OPTEGO / USA
3391 E Hemisphere Loop
Tucson, AZ 85706
Phone 866-678-6723
Fax 520-209-1586
email@optego.com
www.optegousa.com

OPTEGO / CANADA
3600B Laird Rd, Unit 3
Mississauga ON L5L 6A7
Phone 1 416-780-1289
Toll Free 1 877-678-3464
Fax 1 877-471-6001
orders@optego.com
www.optego.com

OPTEGO / SINGAPORE
60 Paya Lebar Rd #11-33
Paya Lebar Square
Singapore 409051
Phone +65 6742 5940
asia@optego.com
www.optego.com.sg